



WP
Winpact



2017 BROCHURE

Utility Equipment



Winpact Chiller



WCC-100
WCC-101



Cat. No.	WCC-100	WCC-101
Display	LED display	
Power	1000W	
Bath capacity	10L	
Chamber Material	SUS#304	
Controller	PID control, PT100 sensor	
Temperature control range	0°C to +100°C	
Temperature stability	± 0.5°C at 20°C	
Temperature accuracy	± 0.5°C at 20°C	
Rated voltage	110V; 60Hz, 14.5A	220V; 50Hz, 8A
Safety device	<ul style="list-style-type: none"> • Self-diagnosed abnormality display • Electronic overheating thermal fuse protection in increments of 0.1°C • Delayed resume compressor protection 	
Circulation volume	Cycles up to 5.5L/min cooling system with delayed resume protection after power outage	
Dimension (cm) (WxLxH)	Approx. 34x56x67	

Laboratory Bioreactor / Fermentor

Winpact Parallel (FS-05 Series)



Winpact One (FS-06 Series)



1L Double Jacketed Vessel

FS-05

- Control up to eight (8) systems from a single interface



10L Single Wall Vessel
with Heating base unit



1L Double
Jacketed Vessel



1L Single
Wall Vessel



FS-06



Winpact Evo (FS-07 Series)

- Control up to sixteen (16) systems from a single interface



FS-07



5L Air Lifter Vessel



5L Single Wall Vessel
with Heating blanket



Front view



Side view

*0° - 90° rotation,
120° for harvest

- Fully integrated system specifically designed for solid-state fermentation research involving saccharification, hydrolysis and more
- Programmable angle-adjustable vessel tilting and stirring mechanism permits superior sample homogeneity
- Reduced stickiness design of impellers ensure sample integrity throughout the fermentation run
- Integrated motor shaft & air sparger assembly creates precise, disturbance-free controls of aeration and mixing
- Chemically resistant double jacketed borosilicate glass vessel design
- Customizable impeller and aeration controller available

*For more information, please contact your local distributors.

US Office

19959 Sea Gull Way
Saratoga, CA 95070
U.S.A.
T/ +1-408-366-9866
F/ +1-408-446-1107

Taiwan Office



No. 37, Wuquan 5th Road,
Wugu Dist., New Taipei City 24888
Taiwan
T/ +886-2-2298-1055
F/ +886-2-2299-7871

Shanghai Office

Room 612, International business exhibition center,
9300 Hunan Road, Pudong, Shanghai, China
National toll-free No.:400-823-9177
T/ +86-21-50795277
F/ +86-21-50795277

India Office

No.12-13-104, St. No.3, Lane No.1,
Tarnaka, Secunderabad - 500017
Telengana, India
T/ +91-40-27001515
F/ +91-40-27001515

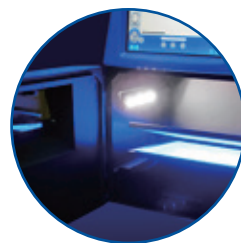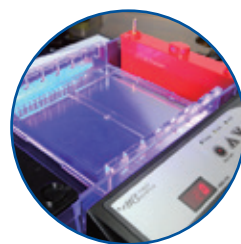

www.majorsci.com
info@majorsci.com